



VSC[®]9000

Redefines the limits of optical
performance to detect fakes,
frauds, & counterfeit documents

foster+freeman



VSC® 9000

Introducing the World's Most Advanced VSC for Document Examination.

Ensuring the authenticity of critical documents is an increasingly complex challenge

From the controlled environment of a police laboratory to the fast-moving world of international border security, the VSC9000 provides a comprehensive forensic-level solution to the interrogation and analysis of *all* document types.

Equipped with a high-performance optical and illumination system, the VSC9000 enables the examiner to conduct multi-spectral examinations of documents in microscopic detail to authenticate legitimate documents, reveal sophisticated counterfeits and to identify evidence of tampering or manipulation with unparalleled accuracy.

Because every microscopic detail matters

The VSC9000 provides a complete forensic toolset to dissect questioned documents, for the inspection of every fine detail from paper and ink composition to the presence of 3rd level embedded security features.

Beyond the physical, the VSC9000 can reveal and analyse every layer of digitally embedded data (including Biometric e-Chip, MRZ, IPI etc.) to provide a comprehensive multi-layered interrogation of a document.



PRECISION
IMAGE CAPTURE



ACCURATE
ILLUMINATION



FORENSIC-LEVEL
ANALYSIS



DIGITAL
AUTHENTICATION

Choose VSC9000 for Forensic-Level Analysis of All Documents.

An evolutionary leap forward in questioned document examination (QDE), the VSC9000 pushes the limits of optical performance for the detection of fake, fraudulent and counterfeit documents.

Choose VSC9000 for Exceptional Image Clarity from Low to High Magnification

Leverages four decades of foster+freeman expertise to capture examination-quality images of documents at up to 480x magnification.

Choose VSC9000 for Unmatched Accuracy in Document Examination

Provides a comprehensive suite of tools for uncompromising detection of counterfeits, forgeries, and signs of tampering or manipulation.

Choose VSC9000 for the Identification of Fakes, Frauds, and Counterfeits

Digitally dissect all forms of document, using advanced forensic analysis tools to unlock hidden evidence, authenticate genuine documents, and reveal fakes and alterations.



TRAVEL
DOCUMENTS




IDENTITY
DOCUMENTS



BANKNOTES
& CURRENCY



OTHER SECURE
DOCUMENTS

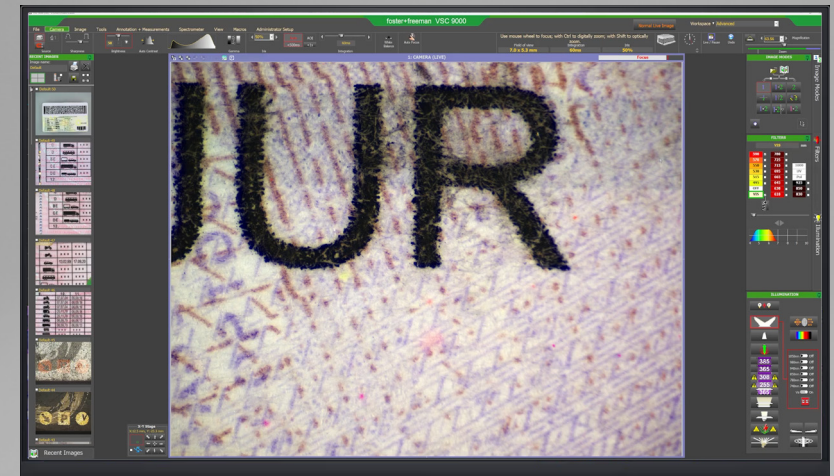

VSC® 9000

Choose VSC9000 to Reveal the Truth Within Every Document

Leveraging over 45-years experience in forensic document examination, the VSC9000 utilises scientific-grade components and is precision manufactured to an exceptional standard to ensure the quality and performance of every system.

VSC Technology

SRI IMAGING

An optical system providing excellent clarity and image resolution for pixel-sharp magnified views of documents.

WaveTRIM ILLUMINATION

Scientific-grade filters refine LED output and provide controlled illumination to reduce glare and enhance contrast and detail for consistent document comparison.

AI Assist MRZ

A machine learning model trained on over 20,000 document images for the automated recognition of document Machine Readable Zones.



12MP
FINE DETAIL



100+
ILLUMINATION
MODES



MICRO
SPECTROMETER



LEVEL 1,2,3
EXAMINATIONS



12MP
OVERVIEW



UV-VIS-IR
LIGHT SOURCES



3D IMAGE
ANALYSIS



E-PASSPORT
ACCESS



480X
MAGNIFICATION



PRECISE
WAVEBAND
ACCURACY



INK&PAPER
ANALYSIS



DATA
DECRYPTION

DUAL-CAMERA
IMAGING

MULTI-SPECTRAL
ILLUMINATION

FORENSIC MATERIAL
ANALYSIS

ADVANCED DIGITAL
INSPECTION

PATENTED TECHNOLOGY

Engineered with cutting-edge technology at our UK facility, the VSC9000 workstations leverage unique patented innovations to deliver unparalleled forensic capabilities.

PAT ID. GB2602793A



UNRIVALLED EXAMINATIONS

A result of over 45 years of VSC design, innovation and refinement, the VSC9000 provides examiners with a comprehensive toolkit specifically designed to address longstanding *and* emerging document security challenges.




DUAL-CAMERA ADVANTAGE

Dedicated optical systems offer seamless magnification transitions between intricate details and full A4/Letter document overviews, maximising efficiency during examinations.



480x IMAGE MAGNIFICATION

Examine in microscopic detail to analyse print and manufacturing quality, and visualise micro-scale security features like taggants.

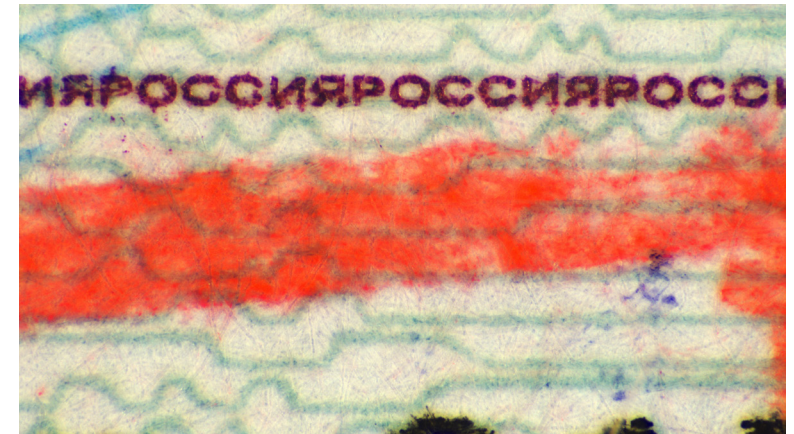

MULTI-LAYERED ANALYSIS

Explore the 3D topography of a document's surface to reveal imperfections, security features, and determine the order that inks, stamps, marks were applied.

MULTI-MODAL ILLUMINATION


Utilise over 100-modes of LED and incandescent lighting to conduct examinations that span the entire forensic imaging spectrum from UV to IR.



VSC[®]9000 Forensic Laboratory

Choose VSC9000 for Legal and Criminal Investigations

For the highest level of document examination, the VSC9000 equips forensic QDE professionals with a comprehensive suite of tools for the dissection and interrogation of all forms of written and printed documents.

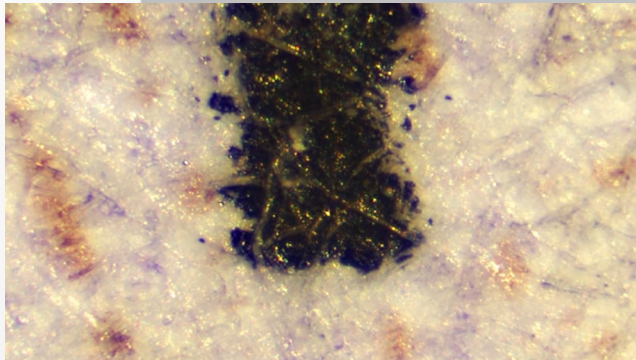


MICROSCOPIC FINE-DETAIL INSPECTION

Reveal alterations and imperfections with pixel-perfect imaging of documents at up to 480x magnification.

Boasting the single greatest leap-forward in image quality during our 40+ year history of VSC product development, the VSC9000 displays live images of unparalleled quality.

The micrometer-precise motorized XY-stage allows examiners to navigate documents and meticulously analyse every detail in real-time.

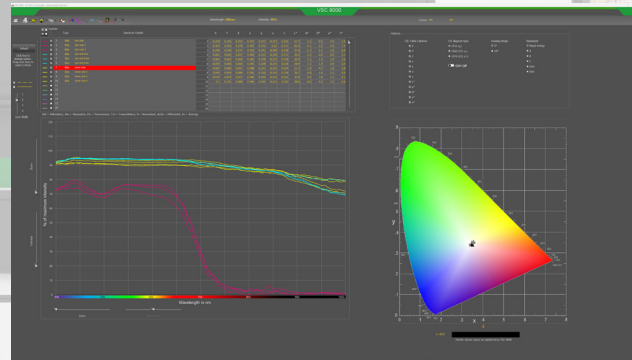

SURFACE & MATERIAL ANALYSIS

Apply scientific techniques to obtain information about a document's material properties.

The system's integrated high-resolution spectrometer delivers non-destructive, quantifiable analysis of a document's material composition.

Measuring absorption, reflectance, transmission, and fluorescence spectra, aids in material identification, detection of alterations, and may support age or origin verification in specific cases.



MULTI-SPECTRAL UV-VIS-IR IMAGING


Exploit light interactions, from UV to IR to examine beyond the limitations of human vision.

The cornerstone of VSC document analysis, precise narrowband illumination can reveal the presence of alterations and enhance the visualisation of security features.

Equipped to examine *any* document, the VSC9000 provides over 100 modes of illumination through a combination of LED and incandescent lamps.

Choose VSC9000 for International Border Security Applications

Providing the final word in document examination for identity authentication, the VSC9000 is an essential tool for border security at air, land and sea ports of entry, immigration centers, and any other location requiring high-level and reliable document checks.

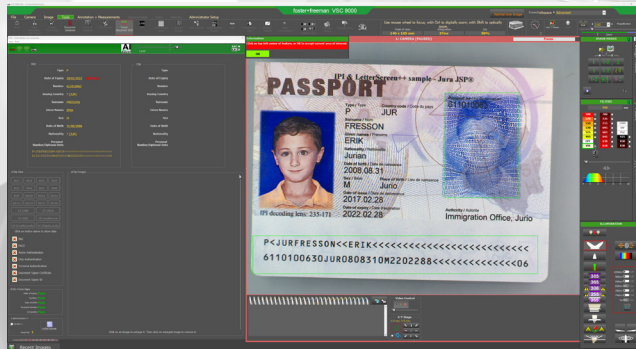


ePASSPORT CHIP DATA ANALYSIS

Extract, decode, and interrogate secure ePassport data using VSC Secure Trust Authentication Chain.

The VSC9000 enables the user to view all personal data contained on the RFID chip of an ICAO-compliant Electronic Machine-Readable Travel Document.


This includes the MRZ and face image protected by BAC and PACE protocols, as well as other biometric and document-specific features.

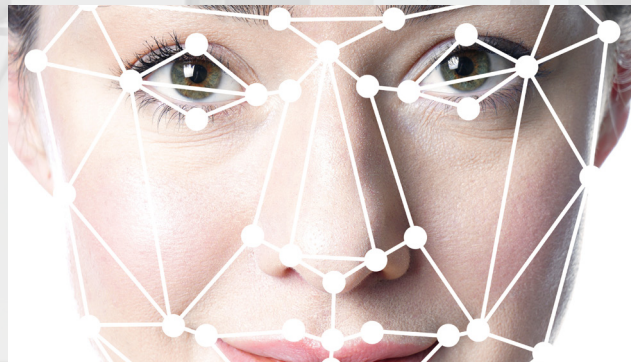

BIOMETRIC FACIAL IDENTIFICATION

Compare and identify document owners using class-leading facial identification technology powered by Thales Multibiometrics.

Equipped with an optional camera, the VSC9000 can provide near instantaneous comparison of live (webcam), embedded (echip), and printed photographs to authenticate the likeness of a document holder. Using the Multibiometrics algorithm ensures efficiency and accuracy for face detection and recognition.

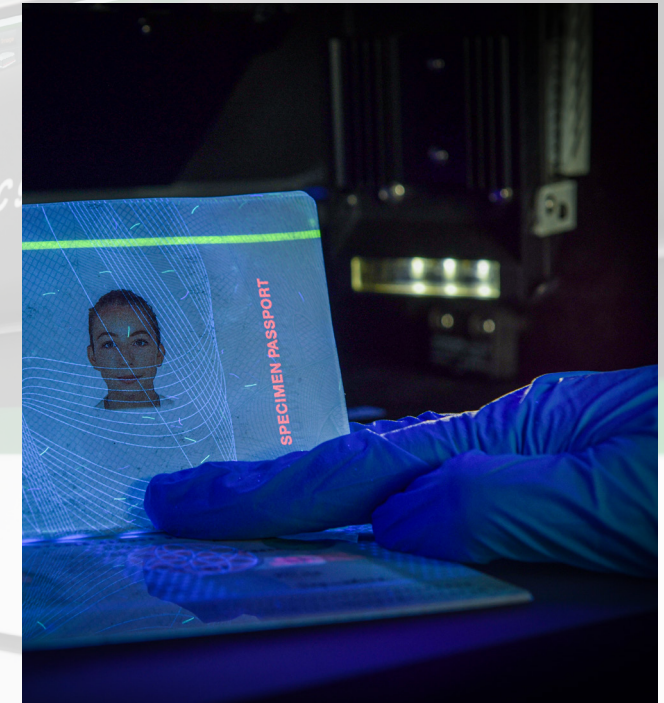

DOCUMENT SECURITY FEATURE ANALYSIS

Activate, visualise and authenticate all levels of security feature printed and/or embedded within a document.

Using a combination of multi-spectral illumination and embedded decoders, the VSC9000 can reveal the the presence of security features including:

Overt | Watermarks, intaglio, security threads etc...
Covert | UV features, complex OVDs, microprint etc...
Forensic | Taggants, nanoprint, phosphorescence etc...

VSC® 9000



HARDWARE features

Optimised for laboratory examinations, the VSC9000 provides imaging, illumination, and material analysis tools for the comprehensive forensic-level analysis of all documents.

Essential VSC Technology includes:

- Dual-camera document imaging system.
- Forensic-level imaging optics providing up to 480x magnification.
- Full-spectral illumination from UV to IR.
- High-specification integral microspectrometer.
- Motorised XY stage for precise micrometre-accurate inspection.

SOFTWARE features

The VSC9000 software provides complete control of examinations and will be instantly accessible to anyone familiar with the VSC line of instruments.

Essential VSC Software functions include:

- Comprehensive document examination toolset with casework management.
- Reveal, decode and display physical and digital security features.
- Extract and examine ePassport biometric personal data.
- Compare documents with stored or database reference documents.
- Advanced analysis modes including 3D and HyperSpectral imaging.

foster+freeman Technology

SRI

SUPER-RESOLUTION IMAGING (SRi)

.. an optical system providing excellent clarity and image resolution of documents captured using a precision 12MP camera.

SRi mathematically processes and combines captured images to create a final Super Resolution Image, which allows for pixel-sharp magnified views of documents. By increasing resolution and clarity, SRi imaging can potentially reveal subtle details in documents that might not be visible with standard magnification.



SEGMENTED COAXIAL LIGHTING ARRAY

Patent GB2602793A

.. a coaxial array providing consistent and controlled illumination along the same axis as the camera lens.

A large, segmented coaxial illumination array allows larger coaxial features to be imaged. The single on-axis segment features improved coaxial contrast, whilst off-axis segments allow sub-surface laser-engraved features to be visualised.



WAVETRIM COLOUR-ACCURATE ILLUMINATION

.. a specialist filter technology that further refines an LEDs output to ensure a high-degree of light purity and waveband accuracy.

Scientific-grade filters 'trim' LED output and provide controlled illumination to reduce glare, boost contrast and enhance detail for consistent document comparison. A unique dichoric coating is applied to each WaveTrim filter to prevent auto-fluorescence and further reinforce the output waveband.



VSC SUITE: AI-ASSIST MRZ

.. for the automated recognition of document Machine Readable Zones (MRZ).

A machine learning model trained on over 20,000 document images, this simple AI tool assists the examiner by automatically recognising the presence of MRZ data within a document. Once detected the MRZ can be read and analysed using Optical Character Recognition (OCR).



VSC SUITE: VSC STAC

.. an optional software module providing Secure Trust Authentication Chain (STAC) analysis of ePassports.

Unique to foster+freeman, VSC STAC enables the user to view all personal data contained on the RFID chip of an Electronic Machine-Readable Travel Document. This includes the MRZ and face image protected by BAC and PACE protocols, as well as other biometric and document-specific features.



VSC SUITE: FF-3D

.. a dedicated software solution for the three-dimensional examination of documents.

Derive a 3D model of a documents surface topography for the examination of intersecting lines, indentations, surface defects, and 3D print features (intaglio etc.). Using FF-3D, visible and IR images can be colour-mapped, inverted, and rotated throughout 3-dimensional space to aid analysis.

VSC® 9000



A Natural Upgrade, for a Step-Change in Optical Performance

17-years since the release of the groundbreaking VSC6000 workstation, the VSC9000 delivers a step-change in optical performance, capable of authenticating a wider range of security features and identifying the high-quality counterfeits currently available via the criminal black market.

	VSC6000	VSC8000	VSC9000
Year of Release	2007	2015	2024
Imaging System	Single-Camera	Single-Camera	Dual-Camera
Image Resolution	5MP	5MP	12MP/12MP
Std. Field of View	210 x 160mm	209 x 167mm	310 x 232mm
Magnification	Up to 170x	Up to 175x	Up to 480x
Microspectrometer	9nm resolution	4nm resolution	2nm resolution
XY Stage	N/A	Motorised	Motorised

VSC9000 System Specification

For full instrument specifications, please contact sales@fosterfreeman.com

Document Imaging System

Primary Camera	- Colour/IR high resolution, live 12MP - Spectral range, 350-1100nm
Secondary Camera	- A4/Letter overview image, live 12MP - Spectral range, 350-1100nm

Standard Illumination Sources

Incident Flood	Visible and IR	LED
	UV-A, UV-B, UV-C	LED
Transmitted	Flood: UV, Visible, and IR Spot: 400-1000nm	LED
Oblique/Side	Visible and IR	LED
Variable Bandpass Spot	- 16x Low Pass filters - 16x High Pass filters	250W
Narrow Bandpass Spot	400-1000nm (1nm steps)	100W

Specialist Illumination Sources

OVD	Two-axis Visible array	LED
Coaxial	Segmented array	LED
Anti-Stokes	IR 980nm	LED
3D Imaging	Visible and IR	LED
Phosphorescence	Pulsed UV	LED
Polarised	Visible Polarised	LED

Microspectrometer

Absorption, reflectance, transmission and fluorescence spectra	Spectral Range: 400-1000nm Spatial Resolution: 2nm
--	---

XY-Stage

Motorised stage	Movement range of at least 60 x 45mm Step resolution 0.01mm (10µm)
-----------------	---

VSC9000 Examination Accessories

For full accessory specifications, please contact sales@fosterfreeman.com

System Accessories

Part ID#

EXTERNAL 32" MONITOR FOR THE VSC 9000 PC	VSC9000/MON
OPTIONAL WIRELESS KEYBOARD AND MOUSE	VSC9000/WKM
OPTIONAL WIRELESS PRECISION MOUSE	VSC/SPACEMOUSE

Examination Accessories

Part ID#

ATR THERMOCHROMIC INK VIEWER	VSC/ATR
POLARISAFE UV LIGHTING MODULE	VSC/POLARISAFE
DOCUMENT TILT PLATE	VSC/TILT
MFSI-2 MAGNETIC SECURITY FEATURE VIEWER	VSC/MSFI/2
QUARTZ GLASS PLATE	VSC/QPLATE
LARGE A4 QUARTZ GLASS PLATE	VSC/QPLATEA4
CARD EDGE VIEWING MIRROR	VSC/EDGEVIEWER
3D IMAGING SOFTWARE MODULE	VSC/3D
THICK ARTEFACT EXAMINATION TRAY	VSC9000/TRANSTRAY

Facial Identification Accessories

Part ID#

FACIAL ID EXPERT HARDWARE KIT	VSC/FR/KIT/E
FACIAL ID STANDARD HARDWARE KIT	VSC/FR/KIT/S
LIVE FACIAL IDENTIFICATION SOFTWARE	VSC/FR10

e-Reader Accessories

Part ID#

EREADER 2A	VSC/EREADER/2A
EREADER 5	VSC/EREADER5
EREADER 6	VSC/EREADER6

Calibration Accessories

Part ID#

NIST TRACEABLE CALIBRATION STANDARDS	VSC/CALIB
--------------------------------------	-----------

ACCESSORIES VSC® 9000

Reference Database Subscriptions

Part ID#

KEESING TRAVEL AND IDENTITY DOCUMENTS DATABASE (UPDATES SUBSCRIPTION)	VSC/DB/KDATA
KEESING TRAVEL AND IDENTITY DOCUMENTS DATABASE (ARCHIVE)	VSC/DB/ARCHIVE
KEESING BANKNOTES DOCUMENTS DATABASE (UPDATES SUBSCRIPTION)	VSC/DB/KDATA/C
KEESING BANKNOTES DOCUMENTS DATABASE (ARCHIVE)	VSC/DB/ARCHIVE/C
ASSURE ID SENTINEL (5-YEAR SUBSCRIPTION)	VSC/ASSUREID-SENTINEL

Embedded Data Decoders

Part ID#

VSC SECURE TRUST AUTHENTICATION CHAIN (STAC)	VSC/STAC
EXTENDED RFID CAPABILITY (SECUNET SDK)	VSC/RFID
LETTERSCREEN ++ DECODER	VSC/LS/PLUS
IPI DECODERS	VSC9000/IPI

Microscope Accessories

Part ID#

EXTERNAL CAMERA (FOR MICROSCOPE INPUT)	VSC9000/CAM
F+F DVM DIGITAL VIDEO MICROSCOPE	V6/DVM2000
NIKON SMZ 745T MICROSCOPE (BASIC)	VSC/NIKON/SMZ745T
NIKON SMZ 18 MICROSCOPE (INTERMEDIATE)	VSC/NIKON/SMZ18
NIKON SMZ 25 MICROSCOPE (EXPERT)	VSC/NIKON/SMZ25
LEICA IVESTA 3 MICROSCOPE (BASIC)	VSC/LEICA/IVESTA3
LEICA M80 MICROSCOPE (INTERMEDIATE)	VSC/LEICA/M80
LEICA M205C MICROSCOPE (EXPERT)	VSC/LEICA/M205C
F+F PVM PORTABLE VIDEO MICROSCOPE	VSC/PVM

FOR FURTHER INFORMATION VISIT www.fosterfreeman.com/VSC9000/

VSC®9000

Ready to find out more?



fosterfreeman.com/vsc-9000/

foster+freeman UK
info@fosterfreeman.com
+44 (0)1386 768 050

foster+freeman USA
usoffice@fosterfreeman.com
888 445 5048

foster+freeman EU
de.sales@fosterfreeman.com
+49 (0)210 255 76 525

foster+freeman NL
nl.sales@fosterfreeman.com
+31 (0)6 1114 3858

foster+freeman